



**FACE & BODY REMODELING**


THE ONLY AUTO ADJUSTING THERMAL SKIN TREATMENT FOR COLLAGEN REMODELING AND SKIN CONTRACTION FOR THE FACE, NECK AND BODY AREAS.

FORMA™

PLUS™



Two configurations of non-invasive sub dermal heating for collagen remodeling. RF current flows uniformly between the electrodes to provide a comfortable thermal experience. Built in temperature controls allow for the attainment of clinically proven optimal epidermal temperatures, and the exclusive capacity for prolonged exposure at therapeutic temperatures.

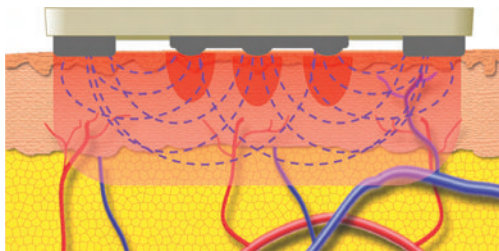

● Heat  
● RF Energy

**A.C.E. TECHNOLOGY BENEFITS**

- **Acquire:** Acquire: Smart RF thermistors built into the hand piece reads skin temperature 1000x per second so clinicians can acquire optimal epidermal settings.
- **Control:** Unprecedented safety with programmable "cut off" temperature.
- **Extend:** Clinical evidence indicates prolonged exposure to (41°C-43°C) is advantageous for skin tightening.

**ROI** Most other manufacturers sell body contouring as a stand alone \$80K+ system –BodyFX is one of your InMode modalities for an all-in-one aesthetic solution... with no consumables!

**SPECIFICATIONS**


RF Output Frequency	Up to 65 W
Output Frequency	1MHz
Temperature Cut-off	Real time temperature measurements set by operator up to 43°C



Dr. H. Ohanian - Treatment of Lax Skin on Skin Type VI



Dr S. Mulholland - After 6 Forma Sessions at 50 mj

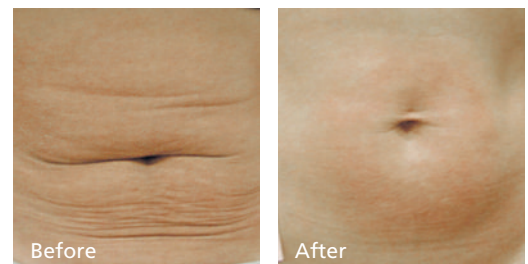

Dr. J. Shaoul - After 6 Treatments at 1 Month Post



Dr. S. Mulholland - 3 Months After 8 Treatments