

# LUMECCA™

by Inmode

## IPL SKIN REJUVENATION

BREAKTHROUGH IPL THAT DELIVERS 3X MORE ENERGY IN THE 500-600NM RANGE TO IMPROVE EFFICACY FOR RED & BROWN SPOTS.




Optimized for light and dark skin, Lumecca increases the potential for clinicians to treat a variety of skin types and conditions ranging from facial pigmentation, superficial vessels, skin texture and/or photodamage.

### BENEFITS

- High efficiency of vascular and pigmented lesion due to high peak power
- Reduction in clinician treatment time from large treatment spot size and high pulse repetition rate
- Improved patient comfort (no pain) and safety from strong sapphire cooling tip
- Fast pulse repetition rate at all settings

### ABSORPTION COEFFICIENT CHART

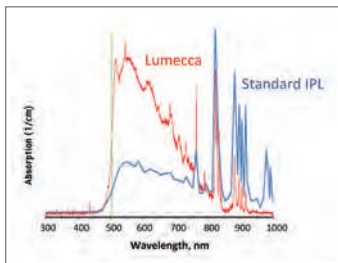

Lumecca treats a wide array of skin rejuvenation concerns with a novel designed low pressure xenon flash lamp with the highest peak power that produces efficient energy at 40% energy in the 500-600 nm range (while standard IPL's only deliver 10-15%).

### SPECIFICATIONS

Wavelength	515 SR	515-1200 nm
Wavelength	580 SR	580-1200 nm
Spot Size		30 mm x 10 mm
Fluence		10-30 J/cm <sup>2</sup>
Pulse Duration		5-30 msec (Short/Long)
Light Cooling Guide		10°-25°C
Repetition Rate		1 pps



Before



After

Dr S. Mulholland - Improvements in brown age spots



Before



After

Dr S. Mulholland - Improvements in brown age spots



Before



After

Dr S. Mulholland - Suitable for multiple treatment areas

